

Measure growth patterns and align sampling spots in photographs: *sclero* tutorial

R package version 0.3.4

Mikko Vihtakari

June 17, 2019

Contents

1	Short introduction to the package	2
1.1	Installation	2
1.2	Features	2
2	Photographs to <i>sclero</i>	3
2.1	Prerequisites	3
2.2	Annotation of photographs in ImageJ	4
2.3	Reading ImageJ zip files containing ROI objects	7
3	Measuring growth	9
4	Aligning sample spots with growth lines	10
4.1	Aligning sample spots	11
4.1.1	<code>spot.dist</code> settings	14
4.1.2	<code>spot.dist</code> troubleshooting	15
4.2	Estimating spatial averaging error of sample spots	16
4.2.1	Add sample spot size information in ImageJ	16
4.2.2	Calculating shell spot size and maximum extent along distance axis	17
4.3	Generating sample maps	19
4.4	Aligning multiple sample spot sequences	21
5	Graphics used in the <i>sclero</i> package	22
6	Dependencies	22

1 Short introduction to the package

The *sclero* package for [R](#) ([R Core Team, 2019](#)) provides tools to measure growth patterns and align sampling spots in chronologically deposited materials. The package is primarily intended for the fields of sclerochronology, geology, and dendrochronology, but the functions can also be applied for other image measuring tasks. The package is developed to interact with [ImageJ](#) ([Schneider et al., 2012](#)), a free-to-use public domain software for image processing and analysis. Since the *sclero* package is developed for open source software, it is free to use and the code is modifiable by anyone interested. If these modifications are distributed in other packages or software, references to the original source (type `citation("sclero")` in [R](#)) and author(s) of the particular function are required. If you use the package in a scientific publication, please cite the package, since it is written as a volunteer contribution. The package might contain errors and users should critically evaluate the results of any function before publishing them. Contributions, code-fixes, and bug reports are welcome and should be committed on [GitHub](#). Persons contributing code to the package will be credited with authorship (*Authors* field in the [DESCRIPTION](#) file).

1.1 Installation

A stable version of the package can be found from [CRAN](#). If you want to use the developmental version (possibly with new features), please install it from GitHub using the [devtools](#) package:

```
devtools::install_github("MikkoVihtakari/sclero")
```

1.2 Features

This tutorial, the functions, and the documentation are subject to changes. Currently *sclero* package contains following features:

1. [Measuring growth](#)
2. [Aligning sample spots with growth lines](#)

2 Photographs to sclero

The *sclero* package uses ImageJ for annotation of features in photographs that are required for further calculations. The package understands four types of ImageJ objects: **straight line**, **segmented line**, **multi-point tool** and **oval selection**. How these objects are interpreted by the package is explained in Figure 1. Click [this link](#) for instructions how to use ImageJ.

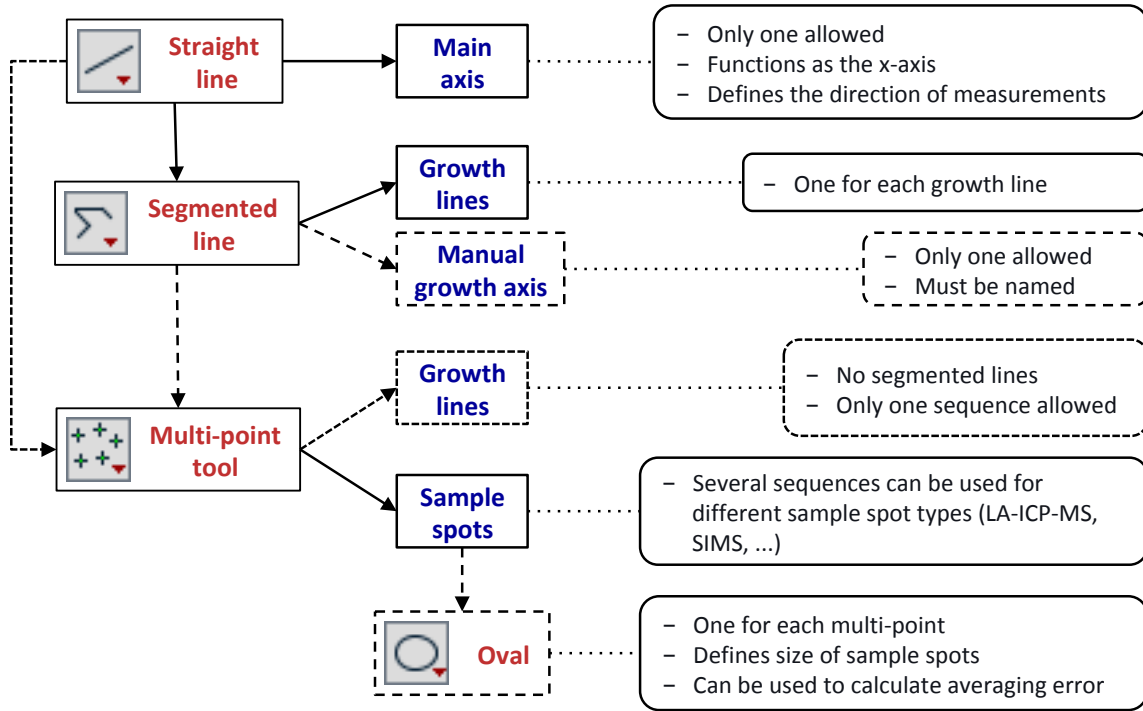


Figure 1: Photograph annotation decision tree. Boxes with red text indicate ImageJ objects within an ROI zip file and boxes with blue text the resulting sclero objects. Boxes with dashed lines are optional. A straight line is always interpreted as the main axis, to which all other objects are scaled. Only one straight line per ROI file is allowed. Segmented lines are generally interpreted as growth lines, but one can be assigned as a manual growth axis. If segmented lines are not provided, multi-points will be interpreted as growth checks for growth calculations. Otherwise, multi-points are interpreted as sample spots. The area of sample spots for averaging error estimation can be defined by oval objects, which must be placed on top of each multi-point.

2.1 Prerequisites

In order to use the *sclero* package, a photograph of sufficient resolution to see the growth lines and other desired features is required (Figure 2). In addition, the user should take care that following points are considered:

Aspect ratio The aspect ratio of pixels in the photograph has to be 1:1, meaning that a certain length vertically equals to the same length horizontally.

File format The image file can be of any format compatible with ImageJ, but .tif images are generally recommended for editing purposes as this format is lossless meaning that it can be saved again without losing information. Yet, as long as the compression does not distort the image, file format should have a minimal impact on the outcome, because the sclero package uses relative distances.

Compiled images If you are using compiled images from a microscope software, pay a special attention to aligning the photographs correctly as this will affect the results. Make sure that all compiled images are taken with a same magnification.

Scaling If you want real distances as output, the microscope has to be calibrated. If real distances are not important, knowing the microscope calibration is not essential as sample alignment is relative to one photograph. Currently there is no way to export the ‘Set Scale’ option from ImageJ and the scaling has to be conducted using the `read.ijdata` function.

Cracks The sample material has to be chronologically deposited without spatial caps or cracks. If your sample material has large caps or cracks, they are bound to affect the output. One option is to edit the photograph in a photo editing program and remove these cracks.

2.2 Annotation of photographs in ImageJ

We use a shell margin of a bivalve shell cross-section from [Vihtakari et al., 2016](#) as an example. The example file can be downloaded from [GitHub in .png format](#). The example shell has one sequence of samples taken with LA-ICP-MS (the large holes, Figure 2) enabling sample spot alignment (see Section 4.4 for multiple sequences). Annotation of growth lines and main axis is enough for growth measurements (Section 3). First, we will mark sample spot centroids using ImageJ, which is not required if you only want to measure growth. Once this is done, we will proceed to annotate growth lines and the main axis. These two annotations are always required.



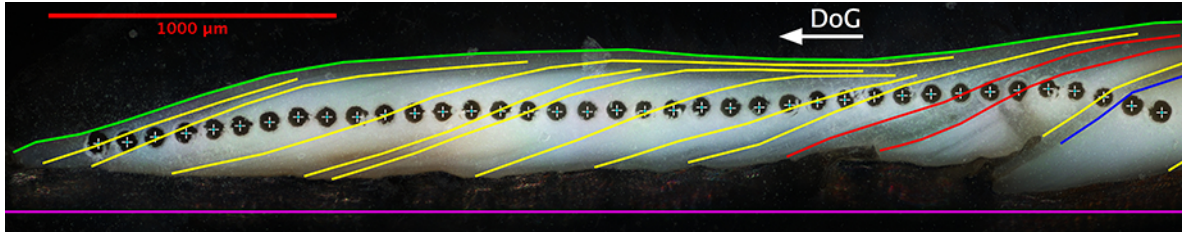


Figure 2: ImageJ annotated bivalve shell section. Colored lines on the shell section indicate the ImageJ marked growth lines, green line further representing the shell margin, yellow lines sub-annual growth lines, red lines the extent of a winter growth band, and blue line a growth line with a known time. Note that the green line (shell margin) covers the entire shell section, whereas other growth lines are long enough to align sample spots chronologically. The magenta straight line represents the distance (main) axis. The main axis should not cross any growth line when `type` is "along". If `type` is "cross", the main axis has to cross all growth lines (see Section 4.1.1). Large holes marked with crosses represent LA-ICP-MS sample spots. White arrow marked 'DoG' points towards the direction of growth. Colors are used here only for illustration purposes and are ignored by the *sclero* package.


1. Start by using the 'Multi-point Tool' and mark LA-ICP-MS sample spot centroids (see Figure 2). The sample spots will be numbered in the order you mark the spots. In this example we begin from the margin and mark the spots against the direction of growth. The location of points marked with ImageJ is considered as the accurate location of each LA-ICP-MS spot, unless the spatial averaging-error estimation procedure is used in addition (explained in Section 4.2).

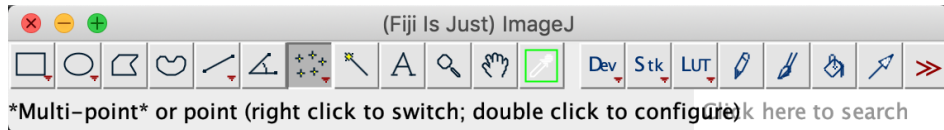

2. After this open the ROI manager and add the sequence using Add[t] button (the keyboard shortcut is `Ctrl + T` or `Cmd + T` depending on the OS). You can rename the ROI (Region Of Interest) to correspond the type of the object ('Laser') for your own reference. The ROI names can be assigned as object names in consequent R functions, although this behavior is not required (see Section 2.3). If you decide to do this, use `_` or `.` as a separator marker. Do not use other special characters (`-`, `,`, `#`) or white space in ROI names, which you want to keep, because these will confuse the internal `grep` functions. ImageJ automatically generates ROI names containing `- (0782-4756` for instance). These names are not valid `data.frame` column names in R and will be regenerated by the functions in the *sclero* package.

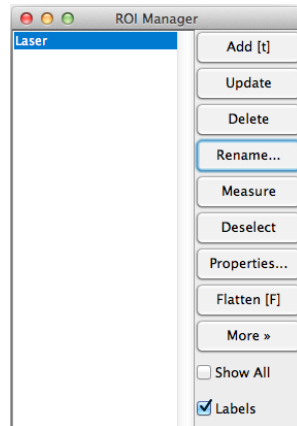

3. After this, mark the visible growth lines using 'Segmented Line' tool starting from the lower margin upwards *against the direction of growth*. Mark each line so that the sample spots are in between at least two lines, but it is not necessary to continue marking much further (Figure 2). Add each line separately to the ROI manager. The order of lines is not important as this will be changed later. Often it is easiest to start with the most clear growth lines and add less clear growth lines in between these lines as needed. Pay special attention that the **growth lines do not cross each other**. Further, there must be **at least one growth line before and after the first and last sample spot** (note the short yellow line on the right in Figure 2). It might be helpful to rename the growth line ROIs, so that they are easier to associate later.

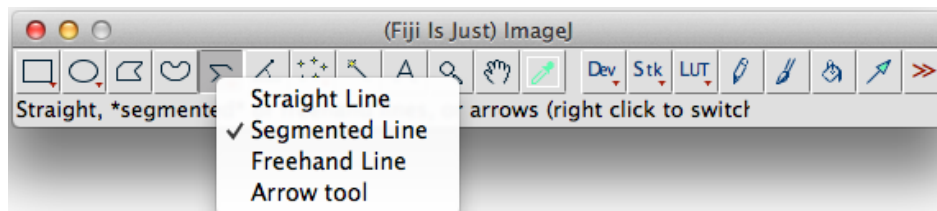

4. Growth measurement functions support one segmented line as manual growth axis. You can add one segmented line as a growth axis if needed. It is a good idea to rename this segmented line as 'growth' for instance.
5. Once you are done with this, add the measurement axis using the 'Straight Line' tool. It is very important to **start the axis from where you want the measurements to begin** as all distances will be scaled to this axis. The main axis should not cross any growth line, if the axis is intended to be approximately parallel to the growth lines (cores and shell sections). If your sample material is a cross-section of an otolith, a tree or a coral or an umbo region of a bivalve, the main axis has to cross all growth lines (see Section 4.1.1). Note that it is possible to have only one Straight Line per ImageJ .zip file since the internal functions recognize the straight line as the main axis automatically.
6. After this you can save the ROI collection as a .zip file (More, Save...) naming it 'shellspots.zip'.

2.3 Reading ImageJ zip files containing ROI objects

Next we can open R and load the *sclero* package.

```
library(sclero)
```

read.ijdata

After this we use `read.ijdata` function to read and process the ImageJ .zip file containing ROIs. The code below is possible to run without having saved ‘shellspots.zip’ file. If you followed the example above, replace ‘path’ with the location of your ImageJ .zip file.

```
path <- file.path(system.file("extdata", package = "sclero"), "shellspots.zip")
dat <- read.ijdata(path, scale = 0.7812, unit = "um")
```

Note that we specified the `scale` and `unit` arguments in `read.ijdata` function, since the example photograph was calibrated. This step is not necessary, but produces actual distances instead of pixel distances (μm in this example; note that directly using μ in R as a character works in some operating systems). The function returns a list of class `IJDATA` containing information about ROIs.

```
summary(dat)
```

	Length	Class	Mode
spots.x	1	data.frame	list
spots.y	1	data.frame	list
gbs.x	14	data.frame	list
gbs.y	14	data.frame	list
main.x	1	data.frame	list
main.y	1	data.frame	list
growth.x	0	-none-	NULL
growth.y	0	-none-	NULL
sample.name	1	-none-	character
scaling.factor	1	-none-	numeric
unit	1	-none-	character
parameters	5	-none-	list

convert.
ijdata

The information in `IJDATA` objects is meant to be user modifiable. The object behaves like any `list` in R only that the name of elements (`$spots.x`, `$spots.y`, ...) should not be changed to ensure that subsequent functions work correctly. In order to align the sample spots along the main axis, we need to convert the information to *spatstat* (Baddeley and Turner, 2005) objects using `convert.ijdata` function.

```
shell <- convert.ijdata(dat)
```

plot.
rawDist

The function returns a list of class `rawDist`, which can be plotted using the generic plotting function. The plot will be a representation of the ROIs, with coordinates rotated relative to the starting and ending point of the main axis.

```
plot(shell, coord.type = "both")
```

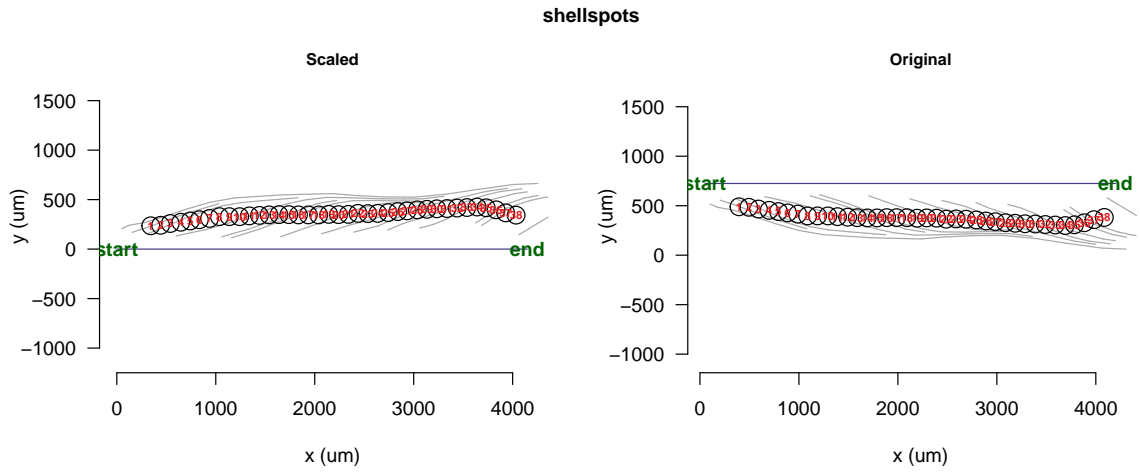


Figure 3: The digitized representation of the shell section. Red numbers show the location of the sample spots. Grey lines represent the marked growth lines and the purple line the main axis to which the sample spots will be aligned. Coordinates are rotated relative to the starting and ending point of the main axis.

3 Measuring growth

Chronologically deposited materials typically contain growth lines, which can be used as time markers to measure growth of the material. Apart from cores, these growth lines rarely are straight or parallel to each other, leading to a situation where placement of the measurement axis (i.e. line or "a ruler" used to measure distance between growth lines) influences the outcome (Figure 4). Typically in such cases, the distance is measured as maximum distance between growth lines (i.e. line of maximum growth; REF).

Figure 4: Bla.

4 Aligning sample spots with growth lines

Materials that can be used as [proxy records](#) in the fields of sclerochronology, dendrochronology and geology are deposited chronologically ([Masson-Delmotte et al., 2013](#)). Consequently, visible growth lines deposited within these records can be used as time markers to reconstruct growth patterns of the material back through time (Figure 5, [Proctor et al., 2000](#); [Sejr et al., 2002](#); [Schöne et al., 2005](#); [Ambrose et al., 2011](#)). In general, these materials do not grow linearly, complicating the dating of geochemical samples from proxy records ([Schöne, 2008](#); [Ambrose et al., 2011](#)). For example, spot samples acquired with laser-ablation inductively-coupled-plasma mass-spectrometry (LA-ICP-MS) or similar techniques along the sampled material are often measured as a distance from a defined position, such as the shell margin in bivalve mollusks, and these distances are then related to time using a variety of geochemical proxies (Figure 5, [Vihtakari, 2014](#)). Sample spots, however, have to be consistently related to the visible growth lines prior time estimates for each sample spot can be acquired (Figure 5). This is because material in proxy records is related to time following the growth rate of the material, and growth lines mark consistent time intervals within these records.

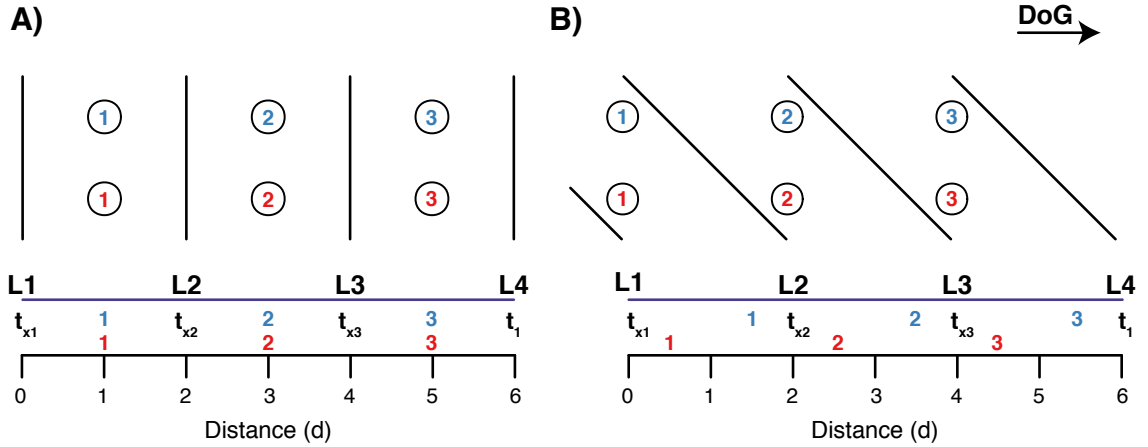


Figure 5: Two sample spot sequences (blue and red numbers inside circles) taken with direction of growth (DoG) along a chronologically deposited material. Growth lines (black lines) mark consistent time intervals (t) relative to a distance from a point with a known time (d) and growth rate of the material (see text). **A)** Growth lines are vertical and both sample spot sequences are perpendicular to the growth lines. Consequently, the vertical position of sample spots does not affect the distance between sample spots and growth lines, and therefore samples taken vertically on top of each other represent the same time intervals (colored numbers without circles above the distance axis). **B)** Growth lines are at a 45° angle to direction of growth, and therefore sample spot sequences are aligned in an angle to the growth lines. Consequently, the sample spots taken on top of each other do not represent the same time intervals. The *sclero* package helps aligning sample spots in both examples correctly using distances from the margin along a distance axis and relative distances of sample spots from adjacent growth lines (see Figure 6).

The importance of the angle between sample spot sequences and growth lines can be illustrated using two simplified examples: the first example (Figure 5A) represents a proxy record where growth lines are at a right angle to sample spot sequences and the direction of growth, such as umbo sections of bivalve shells, tree cross-sections, otoliths, as well as sediment-, rock-, ice- or tree-cores. The three first growth lines from the left (L1, L2, and L3) are deposited after each other at unknown time intervals ($t_x = t_1 - \frac{d}{\bar{r}_x}$, where d is the distance from L4 to the growth line of interest along the distance axis, i.e. $d_{L4} - d_{Lx}$; and \bar{r}_x is the average growth rate of the material between d_{Lx} and d_{L4}). The last growth line on the right represents the margin or the end of the sample material often associated with a known time ($t = 1$ in this example). Calculating the distance of sample spots along the distance scale, which also represents a relative time scale (see the equation above), is straightforward requiring drawing a straight line from sample spots down to the axis leading to distance values of 1, 3, and 5 for both sequences. Some proxy materials, such as margin sections of bivalve shells and heavily bended sedimentary rocks do not, however, allow positioning the sample spot sequences perpendicularly to growth lines. In the second example the sample spot sequences meet growth lines at a 45° angle (Figure 5B). Now both sample spot sequences would yield distance values of 0, 2, and 4, if distance calculation was done similarly to the first example. However, this would lead to a bias when estimating time for each sample spot as the sample spots are obviously not consistently related to the growth lines. To correct for this bias, the sample spots must be related to adjacent growth lines (see Figure 6). This approach gives distance values of 1.5, 2.5, and 5.5 for the blue sequence and 0.5, 2.5, and 4.5 for the red sequence. These distance estimates correspond better with the chronological nature of proxy materials. The distance axis in both examples represents the historical location of the margin allowing direct growth rate calculations, if the time for each sample spot can be estimated.

The *sclero* package together with ImageJ enables alignment of sample spots with visible growth bands from photographs of the sampled material. Despite the simplified example above, the *sclero* package automatically calculates the location of sample spots no matter how non-linear the growth lines are as long as they do not cross at any point. In addition the package estimates the extent of sample spots along the distance axis (see Section 4.2). This extent can be further used in estimation of time-averaging error (Goodwin et al., 2004; Beelaerts et al., 2008). In order to align the samples, the user has to define a distance axis (called ‘main axis’ in the package) to which all sample spots and growth lines will be aligned (Figures 2 and 6).

4.1 Aligning sample spots

`spot.dist`

An object of class `rawDist` can be readily processed using the `spot.dist` function. First, the function projects the beginning of each growth line to the measurement axis (L_1 and L_2 in Figure 6). Then, the sample spots are aligned along the main axis in relation to the adjacent growth lines on both sides of a sample spot such that $d_1/d_2 = d_{L1}/d_{L2}$.

`plot.`
`spotDist`

The function returns a list of class `spotDist` containing new information of the aligned sample spots and the digitized representation of the shell cross-section, which was already included in the `rawDist` object. Also `spotDist` objects can be plotted using the generic plotting command.

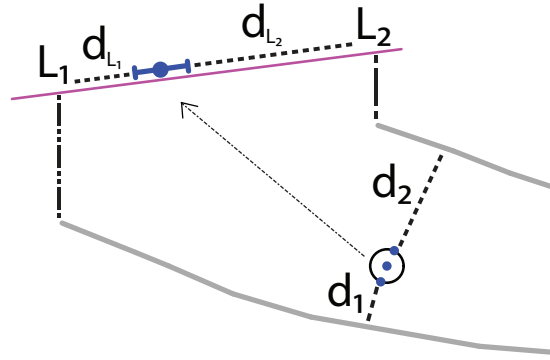


Figure 6: Alignment of sample spots along the distance (main) axis. Grey lines represent ImageJ marked growth lines, open black circle a sample spot, and blue dots the centroid and the closest points to the adjacent growth lines along the perimeter of the sample spot. The `spot.dist` function aligns the blue dots such that $d_1/d_2 = d_{L_1}/d_{L_2}$ resulting to a segment along the distance axis (blue dot with error bars). Only the centroid is used, if sample spot areas are not specified separately (see Section 4.2). The figure is from Vihtakari et al. (2016)

```
# shell <- convert.ijdata(shellspots) # if you did not follow the example above
aligned <- spot.dist(shell)
plot(aligned)
```

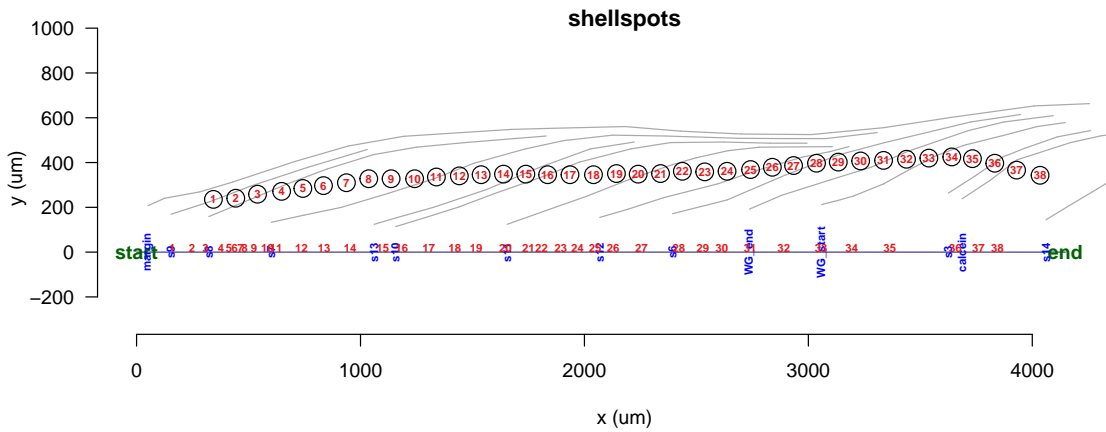


Figure 7: Aligned sample spots along the measurement axis together with the digitized presentation from Figure 3.

Typing the name of a `spotDist` object prints summarized information about the sample spot alignment:

```
aligned
```

```
"spotDist" object for shellspots
```

Main axis type: along

Scaling factor: 0.7812 pixels/um

Including following elements:

	list.elements	type	ncol	nrow	length
main	psp		0	0	5
gbs	psp		0	0	5
spots	list		0	0	1
growth	NULL		0	0	0
window	owin		0	0	4
start.main	ppp		0	0	6
end.main	ppp		0	0	6
sample.name	character		0	0	1
scaling.factor	numeric		0	0	1
unit	character		0	0	1
parameters	list		0	0	10
main.type	character		0	0	1
gb.projections	ppp		0	0	6
gb.start	ppp		0	0	6
gb.end	ppp		0	0	6
gb.projections.dist	data.frame		4	14	4
mid.point.type	character		0	0	1
det.dat	list		0	0	1
output	list		0	0	1

Sample spot output:

spot	dist
1	94.62051
2	182.55956
3	243.07346
4	311.18224
5	346.60488
6	372.90810
7	399.64695
8	417.83668
9	458.90261
10	505.31901
11	544.43780
12	656.68984
13	758.02466
14	874.91661
15	1019.98703
16	1105.41540
17	1226.03512
18	1341.27597
19	1437.40469
20	1567.85498

```

21 1670.48644
22 1729.37731
23 1815.60525
24 1889.99943
25 1968.36939
26 2048.81858
27 2174.97867
28 2341.48465
29 2449.03098
30 2534.59278
31 2660.58426
32 2811.12793
33 2978.10571
34 3115.26282
35 3283.93206
36 3577.47434
37 3680.41824
38 3765.36259

```

Growth line output:

	line	gap	gap.l.main	dist0
margin	margin-s9		102.40655	0.0000
s9	s9-s8		168.97081	102.4066
s8	s8-s7		279.05786	271.3774
s7	s7-s13		459.54941	550.4352
s13	s13-s10		96.00614	1009.9846
s10	s10-s11		497.95187	1105.9908
s11	s11-s12		413.46646	1603.9427
s12	s12-s6		325.14081	2017.4091
s6	s6-WG_end		343.06196	2342.5499
WG_end	WG_end-WG_start		322.58065	2685.6119
WG_start	WG_start-s3		565.79621	3008.1925
s3	s3-calcein		61.44393	3573.9887
calcein	calcein-s14		375.06400	3635.4327
s14	NA		NA	4010.4967

The alignment information with sample spot numbers is stored as a sublist called `output` and can be extracted to a `data.frame` (see below). Otherwise the object behaves like any list in R. Relevant data can be subset as needed. Detailed data containing information of the alignment process is stored in a sublist called `det.dat`.

```

aligned$output ## Results shown above
aligned$det.dat ## Results not shown here to save space

```

4.1.1 spot.dist settings

The alignment function `spot.dist` can either project growth lines on the distance (main) axis as in Figure 6 or use the crossing points between growth lines and the main axis. These

two types of the main axis can be used for different applications. The main axis type is automatically selected by the following criteria:

along Appropriate for samples with cut-off growth lines such as bivalve margin cross-sections and tree, sediment or ice-cores. This option is selected by placing the measurement axis such that **it does not cross any of the marked growth lines**. The location of each growth line is projected along the measurement axis from the beginning of the growth line (the point where you started marking the growth line in ImageJ).

cross Appropriate for approximately round cross-sections: samples where the growth lines continue through the entire width of the sample (such as tree, coral or calcareous algae cross-sections and umbo-regions of bivalves). This type is selected by making the **main axis to cross each individual marked growth line**. The location of each growth line along the main axis is considered as a crossing point.

These criteria are set due to the need of defining a location for each marked growth line along the distance (main) axis. The choice is rigid, to simplify calculations, and to avoid bias in results by allowing two different methods for growth line locations. The easiest way to test which **type** suits a particular sample best is to save two sets of ImageJ zip files by moving the measurement axis.

4.1.2 `spot.dist` troubleshooting

Apart from the common issue of placing the main axis a wrong way around (because where you start drawing the line is considered as the **start** point), problems with `spot.dist` failing to find a location for each sample spot are often connected with the lack of marked growth lines on both sides of each sample spot. If this is the case, try drawing more growth lines so that each sample spot is surrounded by them.


All functions in this package are still at an experimental stage. They do work for the applications needed by the author, but might not work for other applications. They are also likely to contain bugs, which can be fixed. Please contact the package maintainer, if you encounter unexpected behavior or obvious errors in the functions.

4.2 Estimating spatial averaging error of sample spots

The *sclero* package can estimate the size of each sample spot, if size information is available. Size information can be assigned using the ‘Oval’, ‘Elliptical’, and ‘Rectangular selections’ tools in ImageJ. To illustrate the process, we use the example above.

4.2.1 Add sample spot size information in ImageJ


1. Open the *shellspots.png* and *shellspots.zip* files you saved in the previous steps. Both of these files are available in the *extdata* folder of the *sclero* package.
2. Select the ‘Oval selections’ tool and mark the outline of the first sample spot. You can get circles by pressing down *Shift* button while you do the marking. After releasing the mouse, you can still adjust the position of the circle. Once you are satisfied with the size and position of the circle, add the selection to ROI manager (*Ctrl + T* or *Cmd + T* depending on the OS).

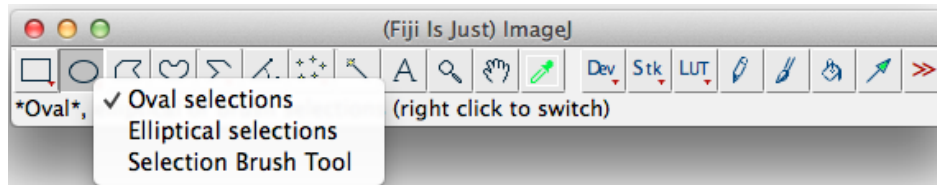

3. Then proceed to the next sample spot and mark it similarly. Mark all the sample spots within the shell section. Note that the order you mark the sample spots does not matter as the *sclero* functions will associate centroid of each sample spot with individual sample spots in a sample spot sequence (*spots* in a *rawDist* object). The ‘Multi-point’ tool marks, however, have to be inside each circle for the routine to work: the perimeters of spot size selections can cross each other, but the centroid of each selection has to be kept adjacent to a ‘Multi-point’ tool mark.



4. Once you have marked all the sample spots using the ‘Oval selections’ tool, save the content of the ROI manager as shellspots.zip. You can overwrite the old file (in fact, if you examine the shellspots.zip file included in the package, you will find out that it already contains spot size information).

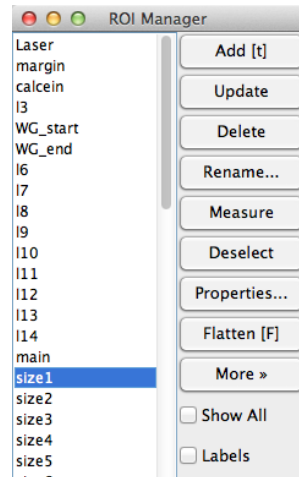

4.2.2 Calculating shell spot size and maximum extent along distance axis

We can either re-read the shellspots.zip file into R (see Section 2.3) or use the `rawDist` object (`shell`), which has already been loaded to R. In here we assume that the user has followed the steps above and has `shell` object loaded in R. If not, run the following code:

```
data(shellspots)
shell <- convert.ijdata(shellspots)
```

Assigning the sample spot sizes to `rawDist` objects is done using the `assign.size` function. If the .zip file containing spot size information is the same than from which the `rawDist` object was derived from and located in your working directory, assignment of spot sizes is simply specified by `assign.size(name_of_the_rawDist_object)`. This tutorial, however, uses the datasets included in the `sclero` package and the path of the shellspots.zip file has to be specified.

```
path <- file.path(system.file("extdata", package = "sclero"))
shellsizes <- assign.size(shell, file = "shellspots.zip", path = path)
```

Sample spot area and diameter can now be extracted from the `rawDist` object.

```
head(shellsizes$spot.area$spot.dat[[1]])
```

Hyperframe:

	spot	dist2spot	spot.owin.name	spot.owins	spot.area	spot.diameter
1	1	4.615401	size1	(owin)	5267.895	83.34566
2	2	3.263581	0001-0378-0383	(owin)	4476.939	76.83428
3	3	3.732048	0001-0365-0459	(owin)	4178.562	74.22973

4	4	1.280161	0001-0355-0546	(owin)	4326.473	75.53200
5	5	1.280003	0001-0343-0620	(owin)	4326.473	75.53200
6	6	1.280082	0001-0334-0693	(owin)	4033.232	72.92745

dist2spot

The `dist2spot` column specifies the shortest distance between recalculated spot centroid location and the sample spot location assigned using the ‘Multi-point’ tool. This distance is used to associate each sample spot circle with a sample spot.

spot.dist
extent

Calculating the maximum extent of each sample spot along the distance (main) axis (Figure 6) can now be conducted using the `spot.dist` function. Centroids of marked spots are automatically used to estimate the location of sample spots along the main axis, and the results therefore slightly differ from those obtained in Section 4.1.

```
spotsizes <- spot.dist(shellsizes)
head(spotsizes$output[[1]])
```

	spot	dist	dist.min	dist.max	spot.area	spot.diameter
1	1	99.74392	22.62012	102.2916	5267.895	83.34566
2	2	180.81162	102.83796	271.3649	4476.939	76.83428
3	3	244.93171	130.33749	271.8315	4178.562	74.22973
4	4	313.45968	270.66711	382.9403	4326.473	75.53200
5	5	344.50228	280.94101	407.8472	4326.473	75.53200
6	6	373.51313	317.17724	429.7806	4033.232	72.92745

The output element of a `spotDist` object now contains not only distance information on sample spot centroids (`dist` column), but also information on the extent (`dist.min` and `dist.max`), as well as spot area and diameter. The actual size of sample spots can now be plotted after a specification.

```
plot(spotsizes, spot.size = "actual")
```

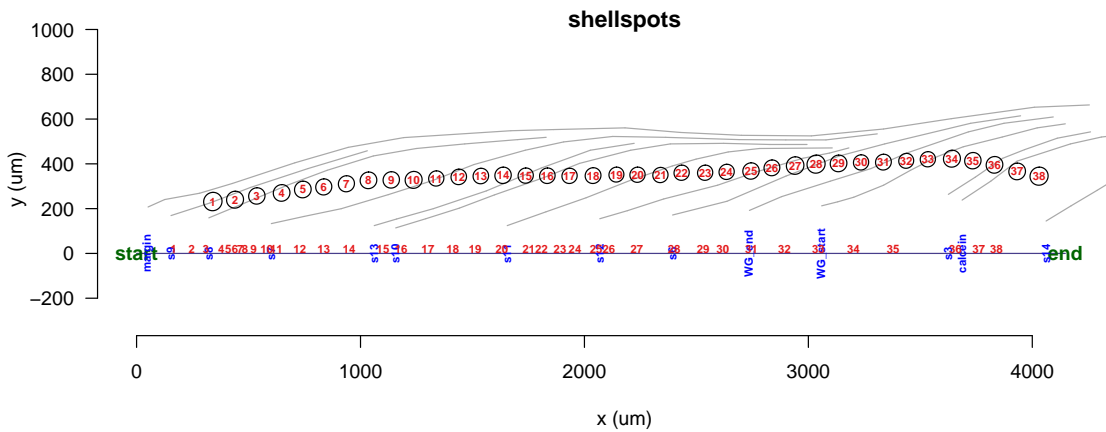


Figure 8: Aligned sample spots along the measurement axis together with the digitized presentation from Figure 3. The circles represent the actual sizes of sample spots.

4.3 Generating sample maps

Sometimes there is a need for plotting sampling results spatially in relation to sample spots. The *sclero* package provides a function for such sample map plotting. As an example we use Ba/Ca ratios calculated from LA-ICP-MS measurements. The dataset is included in the package.

```
data(barium)
head(barium)
```

	spots	Ba
1	1	1.3632722
2	2	1.2361042
3	3	1.1321925
4	4	1.0114307
5	5	0.8961586
6	6	0.7422160

We load the `rawDist` object `shellsizes` from the example in Section 4.2

```
data(shellsizes)
```

`assign.value`

The assignment of values is conducted using the `assign.value` function, where `value` is a data frame with a value for each sample spot. NA's are allowed values and lead to white sample spots. The `spot.type` argument has to be "value" or "idvalue" for density plotting.

```
shellvalues <- assign.value(shellsizes, barium, value.name = "Ba/Ca")
plot(shellvalues, spot.size = "actual", spot.type = "value", main.type = "none")
```

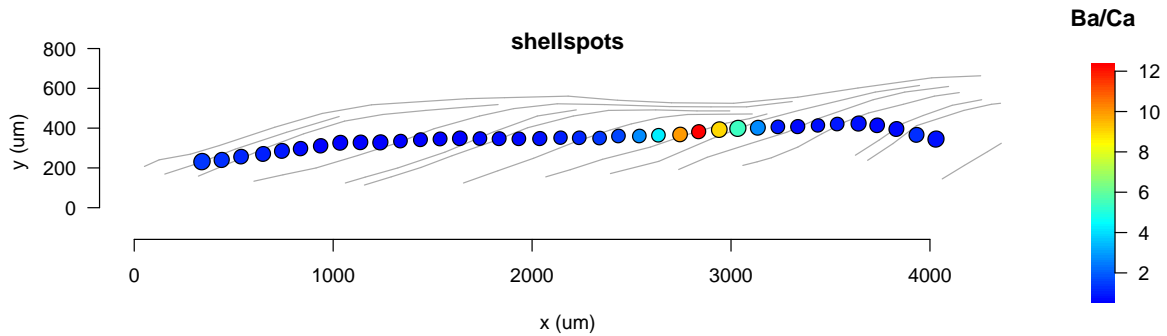


Figure 9: Spatial density map of Ba/Ca over the shell sequence. The `spot.type` argument has to be "value" or "idvalue" to generate such maps. The `main.type` argument is used to remove the distance axis.

A `rawDist` object with value information can be run through the `spot.dist` function and plotted like any `spotDist` object

```

shellvalues.aligned <- spot.dist(shellvalues)
plot(shellvalues.aligned, spot.size = "actual", spot.type = "idvalue",
      spot.color = "darkgrey", highlight.gbs = c("WG_start", "WG_end"))

```


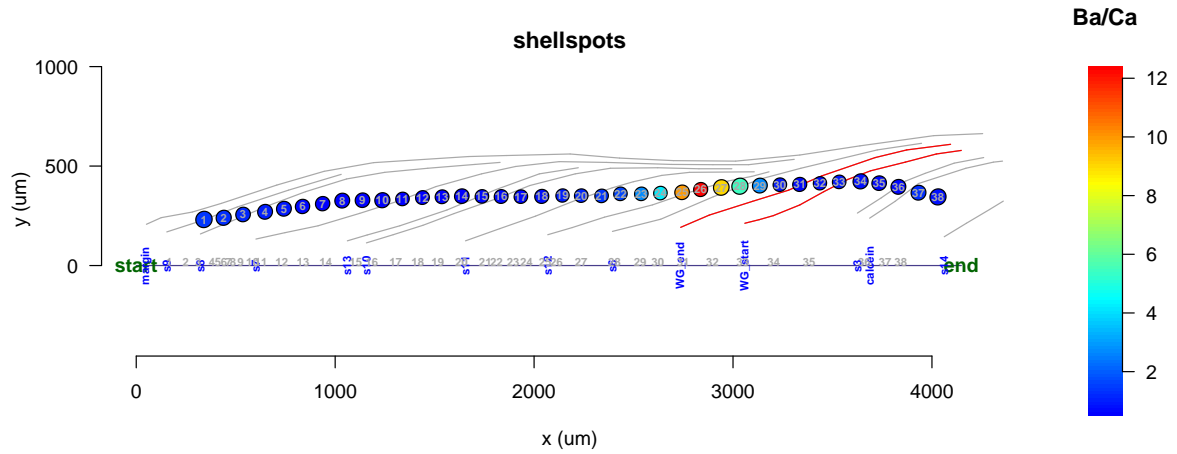


Figure 10: Spatial density map using a `spotDist` object. The beginning and the end of a winter growth band ("*WG_start*" and "*WG_end*") are highlighted and sample spot sequence color specified using the `spot.color` argument.

4.4 Aligning multiple sample spot sequences

Multiple spot sequences can be aligned with the *sclero* package similarly to single spot sequences. These spot sequences must be specified using separate ‘Multi-point’ tool ROI objects for each sequence. In the following example included to the package we use an ImageJ .zip file with a LA-ICP-MS (Laser) sequence and three secondary ion micro-probe (SIMS) sequences. All of the sequences are specified as separate ROIs.

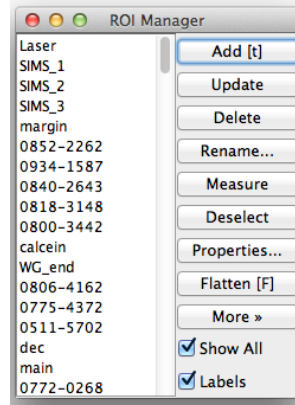

The alignment procedure is similar to single sequences (see Sections 2.3 and 4.1) with the exception that it now produces alignment for each sequence separately.

```
file <- file.path(system.file("extdata", package = "sclero"), "multi_spotseq.zip")
dat <- read.ijdata(file, scale = 0.7812, unit = "um")
multispot.raw <- convert.ijdata(dat)
path <- file.path(system.file("extdata", package = "sclero"))
multispot.size <- assign.size(multispot.raw, path = path)
multispot <- spot.dist(multispot.size)

plot(multispot, spot.size = "actual")
```

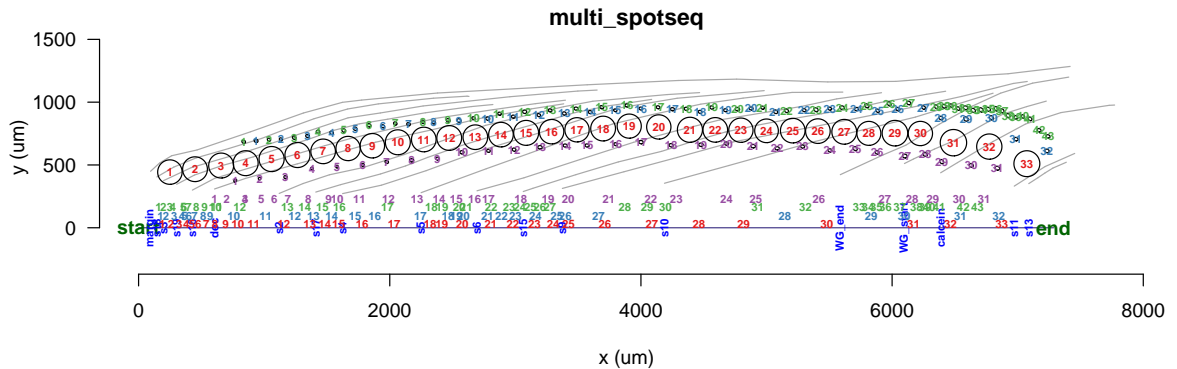


Figure 11: A sample spot sequences with multiple sequences become aligned independently of each other.

5 Graphics used in the *sclero* package

The *sclero* package currently uses the *graphics* package distributed with R for plotting. Plotting sample maps is carried out by the `sclero:::samplemap` function, which works as an internal function and therefore has not been exported. Users willing to modify *sclero* plots beyond the flexibility allowed by `plot.rawDist` and `plot.spotDist` functions are instructed to modify the `samplemap` function, which consists of standard R graphics syntax. It should be noted that `sclero:::samplemap` function calls for the `layout` function every time the arguments `spot.type = "value"` or `spot.type = "idvalue"` are used. Consequently, the graphics window is divided into two regions that might cause issues when combining *sclero* plots with other graphics. The users are advised to consider the graphics window resetting procedure specified in `layout` examples. Any user willing to create a more flexible plotting functionality for *sclero* are asked to contact the package maintainer.

6 Dependencies

The *sclero* package depends on:

RImageJROI (Sterratt and Vihtakari, 2015). Used to import ImageJ ROI objects to R.

spatstat (Baddeley and Turner, 2005). Used for geometric calculations.

plyr (Wickham, 2011). Used for quicker and easier list calculations.

References

- W. G. Ambrose, P. E. Renaud, W. L. Locke, F. R. Cottier, J. Berge, M. L. Carroll, B. Levin, and S. Ryan. Growth line deposition and variability in growth of two circumpolar bivalves (*Serripes groenlandicus*, and *Clinocardium ciliatum*). *Polar Biol*, 35:345–354, 2011. URL <http://link.springer.com/10.1007/s00300-011-1080-4>.
- A. Baddeley and R. Turner. Spatstat: an R package for analyzing spatial point patterns. *J Stat Softw*, 12:1–42, 2005. URL <http://www.jstatsoft.org/v12/i06>.
- V. Beelaerts, F. Ridder, N. Schmitz, M. Bauwens, F. Dehairs, J. Schoukens, and R. Pintelon. On the elimination of bias averaging-errors in proxy records. *Math Geosci*, 41:129–144, 2008. doi: 10.1007/s11004-008-9193-1. URL <http://www.springerlink.com/index/10.1007/s11004-008-9193-1>.
- D. Goodwin, K. Flessa, M. Tellezduarte, D. Dettman, B. R. Schöne, and G. Avilaserano. Detecting time-averaging and spatial mixing using oxygen isotope variation: a case study. *Palaeogeogr Palaeoclimatol Palaeoecol*, 205:1–21, 2004. URL <http://linkinghub.elsevier.com/retrieve/pii/S0031018203007521>.
- V. Masson-Delmotte, M. Schulz, A. Abe-Ouchi, J. Beer, A. Ganopolski, J. G. Rouco, E. Jansen, K. Lambeck, J. Luterbacher, T. Naish, T. Osborn, B. Otto-Bliesner, T. Quinn, R. Ramesh, M. Rojas, X. Shao, and A. Timmermann. Information from paleoclimate archives. In T. Stocker, D. Qin, G.-K. Plattner, M. Tignor, S. Allen, J. Boschung, A. Nauels, Y. Xia, V. Bex, and P. Midgley, editors, *Climate Change 2013: The Physical Science Basis*

- Contribution of Working Group I to the Fifth Assessment Report of the Intergovernmental Panel on Climate Change*, chapter 5, pages 383–464. Cambridge University Press, Cambridge, UK and New York, NY, USA, 2013.
- C. J. Proctor, A. Baker, W. L. Barnes, and M. A. Gilmour. A thousand year speleothem proxy record of North Atlantic climate from Scotland. *Clim Dyn*, 16:815–820, 2000. URL <http://link.springer.com/10.1007/s003820000077>.
- R Core Team. R: A language and environment for statistical computing. R Foundation for Statistical Computing, Vienna, Austria, 2019. URL <http://www.r-project.org>.
- C. A. Schneider, W. S. Rasband, and K. W. Eliceiri. NIH Image to ImageJ: 25 years of image analysis. *Nat Methods*, 9:671–675, 2012. URL <http://www.nature.com/doifinder/10.1038/nmeth.2089>.
- B. R. Schöne. The curse of physiology – challenges and opportunities in the interpretation of geochemical data from mollusk shells. *Geo-Marine Lett*, 28:269–285, 2008. URL <http://link.springer.com/10.1007/s00367-008-0114-6>.
- B. R. Schöne, S. D. Houk, A. D. Freyre Castro, J. Fiebig, W. Oschmann, I. Kröncke, W. Dreyer, and F. Gosselck. Daily growth rates in shells of *Arctica islandica*: Assessing sub-seasonal environmental controls on a long-lived bivalve mollusk. *Palaios*, 20:78–92, 2005. URL <http://palaaios.sepmonline.org/cgi/doi/10.2110/palo.2003.p03-101>.
- M. K. Sejr, K. T. Jensen, and S. Rysgaard. Annual growth bands in the bivalve *Hiatella arctica* validated by a mark-recapture study in NE Greenland. *Polar Biol*, 25:794–796, 2002. URL <http://dx.doi.org/10.1007/s00300-002-0413-8>.
- D. C. Sterratt and M. Vihtakari. RImageJROI: read ImageJ Region of Interest (ROI) files. R package version 0.1.1, 2015. URL <http://cran.r-project.org/package=RImageJROI>.
- M. Vihtakari. *Bivalves as indicators of environmental perturbations related to climate and ocean acidification*. Phd thesis, UiT The Arctic University of Norway, 2014. URL <http://hdl.handle.net/10037/7152>.
- M. Vihtakari, P. E. Renaud, L. J. Clarke, M. J. Whitehouse, H. Hop, M. L. Carroll, and W. G. Ambrose Jr. Decoding the oxygen isotope signal for seasonal growth patterns in Arctic bivalves. *Palaeogeogr Palaeoclimatol Palaeoecol*, 446:263–283, 2016. URL <http://dx.doi.org/10.1016/j.palaeo.2016.01.008>.
- H. Wickham. The Split-Apply-Combine Strategy for Data Analysis. *J Stat Softw*, 40:1–29, 2011. URL <http://www.jstatsoft.org/v40/i01>.